

Bagmatic NOVO




Bagmatic NOVO is our modern blood collection monitor serving the special needs of mobile blood donations. As a light weight device with outstanding long battery life, it is developed to take predefined blood volume and mix the product in three-dimensional way. For quality guarantee, Bagmatic NOVO monitors entire blood donation and registers all significant events that may occur during the process.

Features


- Integrated charger
- Good choice for field applications (battery driven-optional)
- 8 hours of independent work, more than 50 donations autonomy
- 3-dimensional mixing during blood donation
- OLED display with easy menu navigation
- Maximum and minimum flow detection ensures good quality of the product
- Visual indication of error / warning / donation status
- Radio connection for data transfer-optional
- USB cable plug for data collection-standard
- Internal SD card data storage
- Powerful PC software ERYSYS for configuration of the device, data collection and analysis-option
- PC software for collecting and exporting data BNC-standard



Modern intelligent battery concept

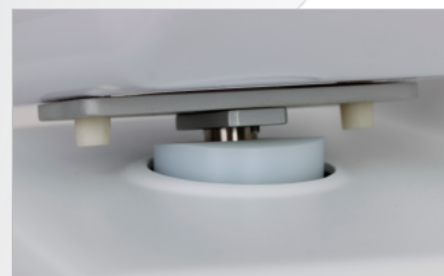

3-dimensional mixing for high quality product

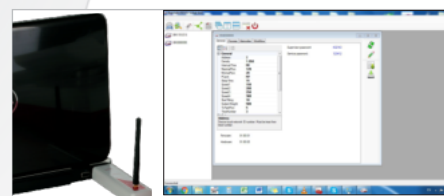

**Simple...User friendly
Real time monitoring
Easy BBIS integration**

Intuitive PC connection software

Data collection


- USB cable connection
- RF Wireless Data communication (option)



Bagmatic NOVO



Illumination ring for visual control of the process

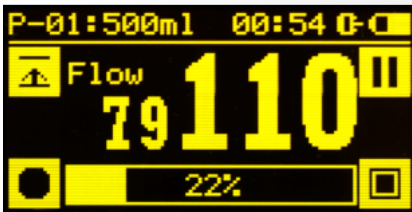

Flow...Volume...
Alarm
Connection status
Battery status,
Time...

OLED display with detailed info

Advantages

- Small size, lightweight
- Designed for mobile long lifetime use
- Battery operated (up to 50 donations or 8 hrs of work) - optional
- Easy to clean and maintain
- Bi-directional communication-option with full software ERYSYS
- Pre-adjusted profiles
- Possibility to install additional language pack
- Customizable alarm system
- Barcode (optional)
- In conjunction with SEALmatic C you can make 2 in 1 device, tube sealer and blood mixer(optional)

Technical Data

Power supply:

AC/DC adapter:	(ACIN: 80-264 VAC; DCOUT: 19 V; 1 A)
Li-Ion Battery:	14,4 V + fast charger (option)
Battery capacity:	3100/3400mAh or more
Battery charging time:	4h

Surrounding temperature:	10-40 °C
Humidity:	maximum 85 %
Device dimensions with Bag tray(L x W x H)	320 x 240 x 215 mm
Weight without battery:	3,55 kg



Bagmatic NOVO with SEALmatic C



**Safe and secure
transport case.
Usable as device stand
(optional)**